

# YPSILON<sup>®</sup> *flow*

*Go with the flow!*

The newest addition to Allard carbon composite AFO technology that is recognised worldwide for providing dynamic response without limiting function.

For active patients who require less stability or more flexibility.





## Features and benefits

- *Optimum range-of-motion in sagittal plane*  
Easily accommodates uneven and sloping terrain as well as longer stride length, which is often a more functional and “normal” gait pattern for active patients.
- *New geometric design for footplate*  
Provides dynamic response for foot lift without restricting or immobilizing normal activity.
- *Lower heel height*  
Increases area of ground contact at stance phase to optimize dynamic response with less effort. Easier for wearer with sagittal plane weakness to “load” the brace.
- *Low-Profile strut connection under footplate*  
Interfaces smoothly with the shoe platform.
- *Gentle forefoot rise*  
The decreased rise in the foot plate provides more space in the toe box, reducing the potential for upper toe box pressure.

## Anatomically-designed tibia interface pads

Evenly disbursed pressure accommodates tibia protuberances. Smooth and thin, this interface pad and strap fits close to the leg, ensuring the brace glides smoothly and is barely noticeable under jeans or trousers. May be applied to accommodate left or right hand pull.



# For your patients that just need a little lift assist!




Ypsilon® Flow is the AFO of choice for your active patients who need sagittal plane assist.

Designed to allow optimum range of motion and dynamic toe-off assist, providing the opportunity for functional or potentially functional muscles, tendons, and ligaments to strengthen and prevent or reduce atrophy.

The footplate is engineered to allow the orthosis to adapt to and move with the lower leg with less resistance to ground reaction forces while the proximal ends of the “Y” provide fixation points, eliminating pressure on the tibia crest.

## Functional Guidelines

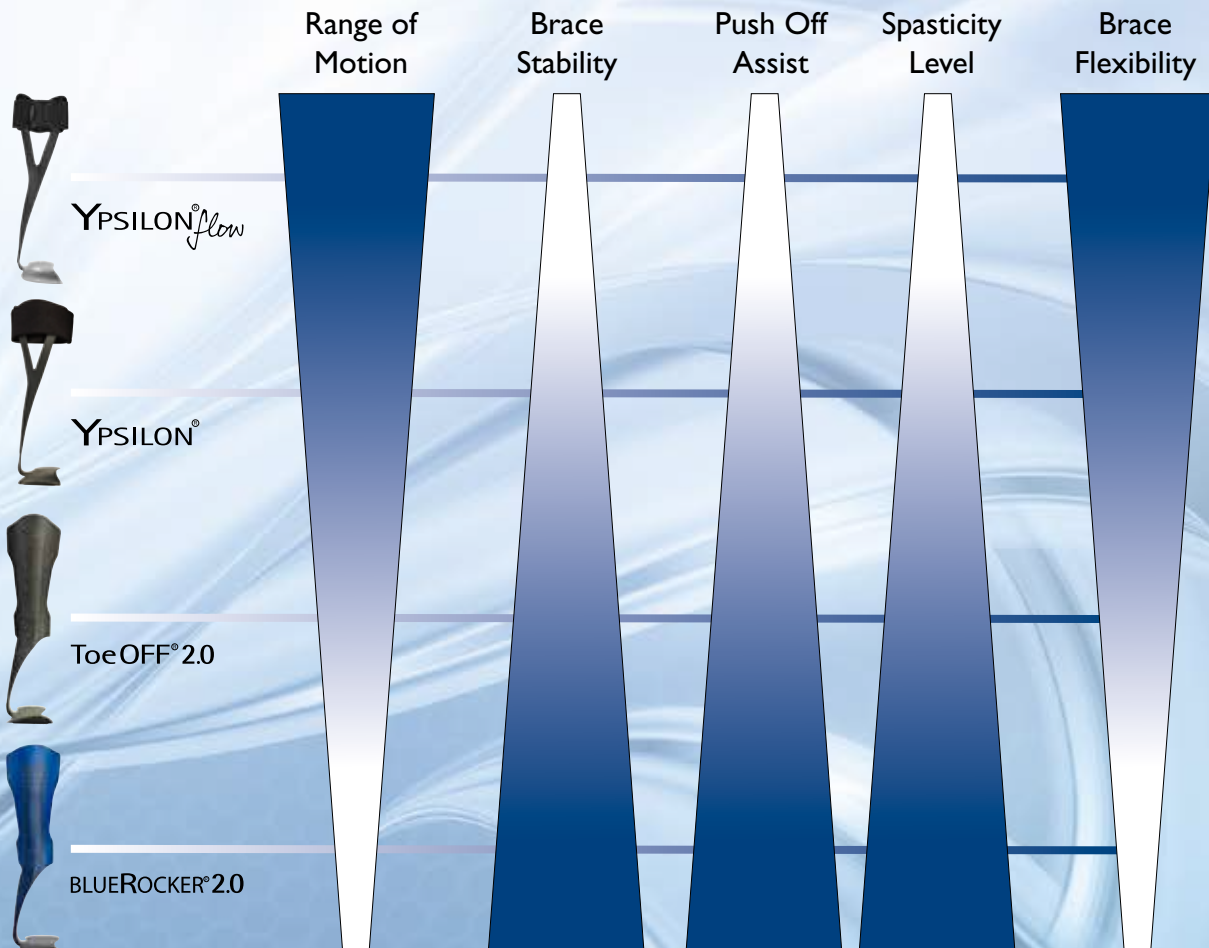

# Support for Better Life!



## Ypsilon®Flow, with Interface and Straps attached

Item no. Left	Item no. Right	Size	Footplate length	Total height	Heel height
28880 1011	28880 2011	Small	230 mm	355 mm	7 mm
28880 1012	28880 2012	Medium	245 mm	365 mm	7 mm
28880 1013	28880 2013	Large	270 mm	375 mm	7 mm
28880 1014	28880 2014	X-Large	285 mm	385 mm	7 mm

## Ypsilon®Flow, Trial Brace and Assessment Kit

Item no. Left	Item no. Right	Size	Footplate length	Total height	Heel height
28881 1011	28881 2011	Small	230 mm	355 mm	7 mm
28881 1012	28881 2012	Medium	245 mm	365 mm	7 mm
28881 1013	28881 2013	Large	270 mm	375 mm	7 mm
28881 1014	28881 2014	X-Large	285 mm	385 mm	7 mm
28882 0000	Gait Assesment Kit, Small - Large, left and right.				

## Ypsilon®Flow, Straps and Tibia Pad

Item no.	Description	Size
28884 0012	Tibia Pad, bilateral fits 3 sizes	Small/Medium/Large
28884 0014	Tibia Pad, bilateral	X-Large
28885 0000	Straps, bilateral fits all sizes	One Size

