


New!



Push button

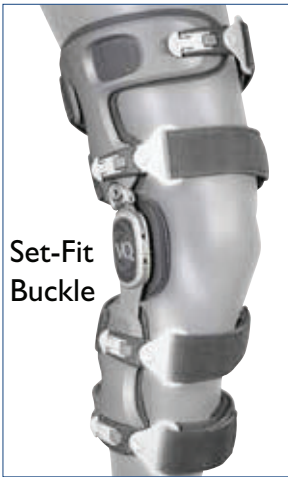

Set-Fit Buckle

OACTIVE®2

OActive® 2 Osteoarthritis Knee Brace is a lightweight, and low-profile upright osteoarthritis knee brace with unique varus / valgus angle adjustment to optimize pain relief. Laterally applied frame and hinge corrects for either medial or lateral compartment and avoids contralateral limb contact during gait. The unique Q-joint is adjusted clockwise or counter clockwise to achieve the desired relief. Start with the orthosis in the neutral position, with the knee bent at 90°, and turn 1/8 turn at a time to adjust the relief. OACTIVE® 2 can also be set to limit the range of motion for extension and flexion independently. The orthosis can also be set for free movement. Comes with 5° stop installed.

Available in Set-Fit or Push Button Buckle System in either:

- Standard hinge for patients weighing under 136Kgs
- Heavy Duty hinge for patients weighing over 136Kgs

Features & Benefits

- Lightweight, low-profile, single upright OA brace
- Unique Set-Fit Buckle System eliminates the hassle of Velcro® adjustments each time the brace is applied and removed.
- Single hand Push Button Buckle for patients with limited hand function
- Maximized 3 points of leverage for medial compartment
- Laterally applied frame corrects for either medial or lateral compartment
- Lateral application avoids opposite limb during gait
- Supplied with Suspension Wrap

Indications:

Mild to severe unicompartmental osteoarthritis.

Item No. Left	Item No. Right	Description	Size
550151000	550152000	Standard Hinge Set-Fit Buckle with Suspension Wrap	<136Kg
550001000	550002000	Standard Hinge Push Button with Suspension Wrap	<136Kg
550161000	550162000	Heavy Duty Hinge Set-Fit Buckle with Suspension Wrap	>136Kg
550171000	550172000	Heavy Duty Hinge Push Button with Suspension Wrap	>136Kg



Suspension Sleeve Under brace



Suspension Sleeve Over brace



New!

TRULIGN OA

Trulign OA™ is a time-tested single upright OActive™ brace design in a patient preferred wrap-style brace. The Q-hinge technology provides significant angular correction to relieve pressure and pain of knee osteoarthritis.

Features and Benefits

Increased patient compliance as it is:

- Easy to apply.
- Breathable, comfortable and conforms to the anatomy.
- Q-hinge angular adjustment delivers optimized 3-point leverage for maximum relief.



Item No. Left	Item No. Right	Size	Circumference 15cm above Patella	Circumference 15cm below Patella
55018 1011	55018 2011	Small / Medium	36-47cm	28-36cm
55018 1012	55018 2012	Medium / Large	47-65cm	36-48cm
55018 1013	55018 2013	Large / X-large	65-81cm	48-64cm

New!

Catalyst Propel OA

This four-point knee orthosis is intended for the treatment of mild to severe unilateral knee osteoarthritis in combination with knee instability or after injury to reduce the load on the joint. The orthosis has the unique Active Thigh Cuff (ATC) solution on the front of the thigh. This is a flexible and strong double thigh cuff that enables the ultimate suspension. Provides protection and comfort by moving the muscles throughout the movement. Strong but low-profile aluminum frame of the highest quality provides superior protection.

The Q-joint of the orthosis is adjusted clockwise or counterclockwise to achieve the desired relief. Start with the orthosis in the neutral position, with the knee bent at 90°, and turn 1/8 turn at a time to adjust the relief. The orthosis comes with different stops for flexion and extension. Catalyst Propel™ OA can also be set to limit the range of motion for extension and flexion independently. The orthosis can also be set for free movement. Comes with 5° stop installed.



Item No. Left	Item No. Right	Model	Size	15 cm below patella	Mid patella	15 cm above patella
55001 1011	55001 2011	Catalyst Propel OA	Small	29-33 cm	29-33cm	29-33cm
55001 1012	55001 2012	Catalyst Propel OA	Medium	33- 37 cm	33-38cm	33-37cm
55001 1013	55001 2013	Catalyst Propel OA	Large	37-42 cm	38-43cm	37-42cm
55001 1014	55001 2014	Catalyst Propel OA	X-Large	42-47 cm	43-48cm	42-47cm
55001 1015	55001 2015	Catalyst Propel OA	XX-Large	47-56 cm	48-56cm	47-56cm